

# YAMAN

TOKYO JAPAN

English User's Manual come with this product.  
Please download the other language's User's Manual version from this link  
(<https://www.ya-man.co.jp/en/support/documents/>) or scan with QR code with  
your camera phone.  
Please carefully read User's Manual to ensure safe and proper use of the device.



TO ACCESS SCAN WITH  
YOUR CAMERA PHONE

## User's Manual

Household Optical Beauty Device  
Rei Beauté Vénus  
Beauty PLUS  
YJEA3

### Introduction

Thank you for purchasing a YA-MAN product.

- Do not use the device before reading the Safety Precautions (pages 3-12).
- Please read this User's Manual carefully for information on how to use the device properly and safely.
- After reading the contents hereof, be sure to keep this User's Manual in an accessible location.
- Make sure to include this User's Manual if you are transferring ownership of the device to another person.

## Steps for Use

1

### Prepare the treatment area

Remove unwanted hairs of the area to be treated with a hand shaver or razor.



2

### Turn on the device

Press and hold the POWER/LEVEL button.



3

### Select the level

Press the POWER/LEVEL button.



4

### Flashing light

Press the FLASH button.



5

### Turn off the device

Press and hold the POWER/LEVEL button.



# Table of Contents

---

Safety Precautions ————— p.3

Names and Functions of Components — p.13

Before Use ————— p.17

How to Use ————— p.20

After Treatment ————— p.28

Additional Information ————— p.30

• Troubleshooting ————— p.30

• Product Specifications ————— p.35

After-sales service ————— p.36

## Safety Precautions

- Read these Safety Precautions carefully before using the device in order to use it correctly.
- Use the device sensibly in accordance with the prescribed steps.
- Immediately discontinue use and arrange to have the device repaired if any abnormality occurs.
- Do not use the device in a defective state, as doing so could cause a physical harm to the user.
- The precautions listed in the User' s Manual are very important. For your safety, please follow them carefully.
- Precautions are classified as DANGER, WARNING, or ATTENTION in order to clearly specify the level of danger and possible injury or damage that may occur if the device is misused.

### Classifications of Precautions



**DANGER**

Improper use of this device poses a severe and likely risk of death or serious injury.



**WARNING**

Improper use of this device may pose a risk of death or serious injury.



**CAUTION**

Improper use of this device may result in an injury or property damage\*.

\*'Property damage' refers to damage affecting a dwelling, household belongings, or domestic animals and pets.

### Symbols



This symbol shows you what you must not do.



This symbol shows you what you must do.



This symbol shows you what you must be careful with.



This symbol shows that strong flashes of visible and invisible light from the device may damage your eyes.

## **WARNING** (The power supply)

- Do not use if the AC adapter is damaged or if the wall socket is too loose.

Failure to comply may cause an electric shock, short circuit, fire, or malfunction of the device.

Contact your local distributor for an advice on repairing the device.



Use of a  
damaged electric  
circuit prohibited

- Always use the dedicated AC adapter.

Failure to comply may cause an electric shock, short circuit, fire, or malfunction of the device.



Prohibition

- Do not exceed the rated capacity of a power outlet or wiring accessory (do not connecting too many plugs in one outlet).

Failure to comply may cause an electric shock, short circuit, fire, or malfunction of the device.



Overloading of  
an electric circuit  
is prohibited

- Do not damage or modify the AC adapter.

Do not apply unnecessary force, bundle, or twist.

Failure to comply may cause the cable to break, resulting in a fire or an electric shock.

Straighten out any twists in the cord before storing.



Careless or  
reckless use  
prohibited

## **DANGER (Use and storage)**

- **The following individuals should not use the device:**

Failure to comply may cause electronic medical devices to malfunction.

◆ Do not use the device in conjunction with any medical electronic equipment, including:

- pacemakers or other implanted electronic medical devices; • artificial heart-lung machines or other electronic life-support medical devices; and
- electrocardiographs or other electronic medical devices worn on the body.



Prohibition

## **WARNING (Use and storage)**

- **Do not use under the following conditions:**

Failure to comply may cause accidents or other problems.

- If you have existing or suspected heart disease or related disorders • If you are unable to sense temperature or stimulation
- Expectant mothers from early pregnancy to just after delivery • If you have an autoimmune disease • If you have a melanoderma
- If you have epilepsy • If you are photosensitive • If you are using photosensitizers, for example if you are undertaking photodynamic therapy
- Infants or those who are unable to express their intentions on their own

Consult a physician prior to use under the following conditions:

- If you have an acute disease • If you have an infectious disease • If you have a malignant tumor • If you have a febrile disease
- If you have high blood pressure or a blood disease • If you have a fever
- If you have sensory impairment due to severe peripheral circulatory dysfunction (e.g. diabetes) • If you have nerve damage
- If you have a skin disease or atopic dermatitis • If you are a hospital outpatient and are taking medication
- If you have a pre-existing skin condition, such as an allergic constitution or contact dermatitis, or if you have a concern about your skin when using this device

- **Do not use in the following cases or on the following areas:**

Failure to comply may cause accidents or other problems.

- If your skin is damaged, experiencing a rash, inflammation, or other abnormalities • If you are under the influence of alcohol
- Anywhere inside the mouth or any mucous membrane • Eyeballs • Eyelids (eye sockets) • Throat • Bulging moles • Viral warts
- Areas where you have broken a bone • Areas of internal bleeding, swelling or scabbing after cosmetic medical procedures
- \* Avoid using the device for 2 weeks after cosmetic medical procedures. After this period, consult with a physician and only use the device if the physician confirms there is no problem. (e.g. hyaluronic acid injections, silicone injections, Botox, thread lifts, double eyelid surgery, HIFU, etc.).
- Dark-colored areas such as scars and spots • Eyebrows, forehead, and hair (including sideburns) • ears, • Nipples, areolae • the nape of the neck
- Inside the navel • Areas where a large vein or other blood vessels can be seen • Tattoo or body-painted area



Prohibition

## **WARNING** (Use and storage)

- Do not overuse the device.

For usage time and frequency, follow the instructions in the User's Manual.

Failure to comply may cause physical injuries.

- Before use, check to make sure that the treatment head, treatment window, and attachment are not cracked, deformed, or broken.

Do not use the device if those parts are damaged or malfunctioning.

Failure to comply may cause injuries or skin problems.



Prohibition

- Do not use, store or charge the device in the following locations:

- Bathroom, shower room, or near the wash basin or other containers (it is okay to use the device)
- A dusty location, a location where there is exposure to direct sunlight or outdoors
- Uneven or unstable location, or a location where there is vibration or shocks
- A location where there are strong radio waves or magnetism
- Near a flammable substance (such as alcohol, thinner, benzene, nail polish remover or spray).
- A location where the device is exposed to oil mist or steam (it is okay to use the device)
- A location that may become hot, such as near a heating device (e.g. a heater)
- A location with drastic temperature changes (for example, near an air conditioner or heater)

Failure to comply may cause a fire, electric shocks, accidents, malfunction of the device, falls or injuries.



Prohibition

## **WARNING** (Use and storage)

- **Do not use the device in hot springs (spas), warm water containing bath salts or sauna.**

Failure to comply may cause a malfunction of the device.

- **Do not leave the device unattended while it is turned on.**

Failure to comply may cause a fire or accident.

- **Do not use near children.**

As strong light is emitted as flashes, they may cause damage to the eyes.

Please be careful not to let people around you see flashes directly.

- **Do not use the device in a dark room. In addition, do not look at, or directly into, the irradiation unit.**

Failure to comply may cause a fire or accident.

- **Do not apply the flash on clothing or objects.**

Failure to comply may cause scorching, discoloration, or damage.



Prohibition

- **Do not direct the treatment window towards a black object.**

Failure to comply may cause a fire or accident.

- **Do not touch the filter (transparent glass) during use and immediately after use.**

Because the filter becomes very hot during use, touching it may cause a burn injury.

- **Do not use the device when the fan is not working (when a humming sound cannot be heard).**

The cooling fan may be broken.

Failure to comply may cause smoke, fire, a malfunction of the device, or an injury.

Please contact your local distributor (page 36).



Prohibition



## **WARNING** (Use and storage)

- Children should obtain the consent of their guardians before using the device and should only use it under the supervision of a guardian or a specialist.

Failure to comply may cause an accident or an injury.



User restrictions

- Keep the device and AC adapter away from children and do not allow them to use it alone.

**Do not allow children to play with the device.**

Failure to comply may cause an accident or an injury.

- Do not disassemble, attempt to repair or modify this device in any way.

Failure to comply may cause a fire, electric shock or other injuries.

Contact your local distributor for an advice on repairing the device.

Do not repair or replace the glass tube.

Do not disable the safety functions (touch sensor etc.).



Disassembly prohibited

- Do not dispose of the device in a fire or place it near fire or flames. Do not use, charge or store in a high-temperature or high-pressure area.

The device features a built-in rechargeable battery.

Failure to comply may lead to the generation of heat, a fire or explosion.



Combustion warning

## **WARNING** (Use and storage)

- When disconnecting the plug from a power outlet, be sure to pull it out by grasping its head instead of the cord (do not touch the metal parts).

Do not insert or disconnect the plug with wet or damp hands.

Failure to comply may cause an electric shock, short circuit, fire, or malfunction of the device.

- If you feel an earthquake, see lightning or hear thunder, turn off the power and unplug the device.

Failure to comply may cause a fire or electric shock.

- After using the device and when cleaning it, be sure to turn off the power and remove the plug from the wall socket.

Failure to comply may cause a fire or accident.

- When not using the product, or when cleaning it, make sure to turn off the power and to unplug it.

Make sure to disconnect the AC adapter from a power outlet after charging is completed.



Removing  
the plug

- Do not leave the device submerged in water or a bath. Failure to comply may cause a malfunction of the device.
- Ensure that the temperature of the usage environment is properly maintained.

The range of environmental temperature in which the device can be used is from 10 to 35°C. Verify the ambient temperature of the location where the device will be used. Failure to comply may cause the device to fail or malfunction.

The device may not correspond effectively to the ambient temperature of a new location immediately after moving to the new location.

- If the device is dropped or damaged, turn off the power immediately, pull out the plug, and request a repair/inspection.

Failure to comply may cause a fire or electric shock.



Must Comply

- Do not use the device for any purpose or in areas other than those stated in the User's Manual.

Failure to comply may cause a malfunction of the device or injuries.



Non-intended  
uses prohibited

## CAUTION (Use and storage)

- Do not allow the device to come into contact with another beauty device or electronic appliance while the power is on. Do not use with other electronic products (electric blankets/bed sheets, heated tables, etc.) at the same time.

Failure to comply may lead to accidents due to malfunctions.

- Do not use while driving.

Failure to comply may cause accidents or injuries.

- Do not press the device strongly or use it intensively on a single area.

Pressing the device stronger than necessary or using it intensively on a single area may damage the skin.

- Do not use the device underwater.

The product is not designed to use underwater, so please be careful when using the device in bathrooms or around water.



Prohibition

- Do not apply flashes on the same area more than once a day

(As for soft flashing, three high-speed continuous flashes are counted as one flash).

- Flashing for an excessively long period of time or in a repetitive manner may cause damage to the skin.

Please observe the guidance in 'Treatment cycle' (page 18).

- There is no problem with the areas treated with flashes overlapping with each other to some extent in each use of the device.

- Do not use perfume, cosmetics, or antiperspirant on the treated areas immediately after treatment is completed (excluding any cream or lotion for moisture retention).

Failure to comply may cause skin problems.



Prohibition

## CAUTION (Use and storage)

- **During the treatment, carefully check that your skin is not showing any abnormalities.**

Stop using the device immediately if you start to feel pain, fatigue or see any skin reaction during use.

Stop using the device immediately if the treatment does not agree with your skin.

If symptoms persist the following day and beyond, consult a physician.

- **Refrain from hard exercise and drinking alcohol on the day of treatment.**

Failure to comply may cause skin problems.

- **For twenty-four hours after the flash treatment, avoid exposing the treated area to ultraviolet rays.**

If redness continues to be exhibited even after 24 hours, avoid UV light.

Also, protect your skin with sunscreen, etc. while undergoing treatment. Failure to comply may cause skin problems.

- **Be sure to wipe off moisture from the main unit after use.**

Also, do not leave the device in hot or humid places, such as in a bathroom.

Leaving the main unit wet could cause damage.

- **Remove any bits of hair remaining on your skin before use.**

Failure to comply may cause your skin to feel hot when you treat the areas with hair debris.



**Must Comply**

## ■ When Using Outside of Japan

### **WARNING**

- Use at AC100-240V only.

Failure to comply may cause a fire or electric shock.



**Must Comply**

- Be sure to check the plug configuration and the voltage of the power outlet before use.

Failure to comply may cause the generation of heat, a fire, or a malfunction of the device.




**Combustion  
warning**

- The plug configuration and voltage may vary depending on the region and the type of accommodation. Always be sure to check for compatibility before using the device.
- For information about plug adapters, consult a consumer electronics retailer.


## Names and Functions of Components

### Main unit

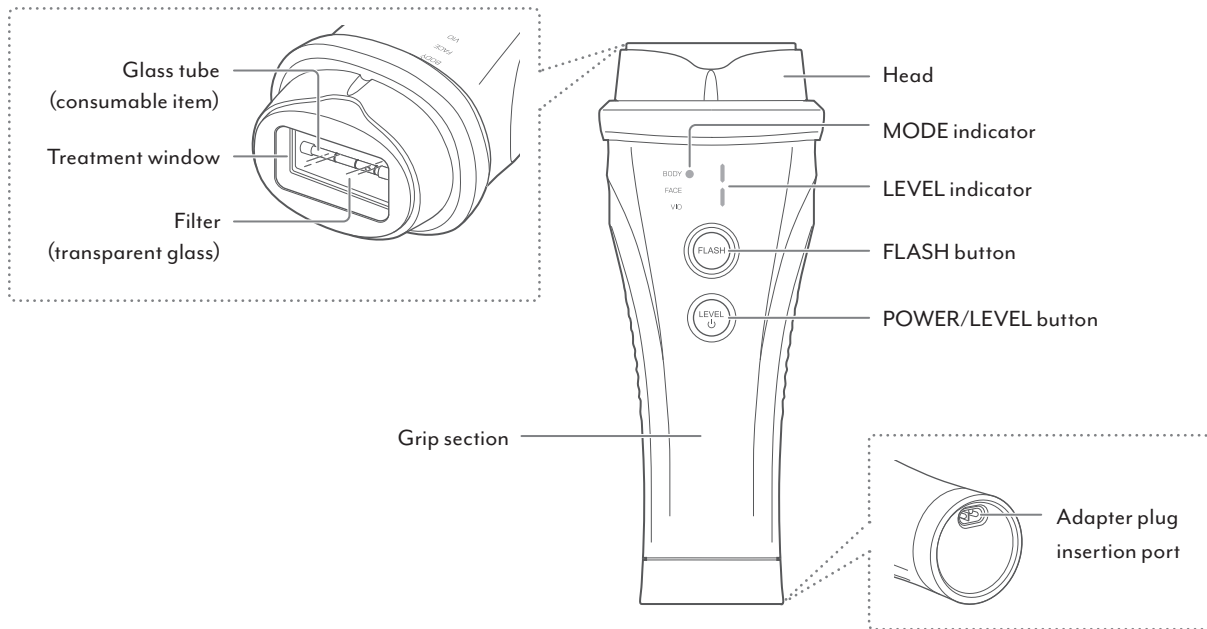


## Accessories

### ● Facial attachment

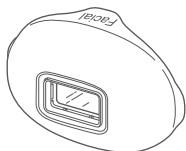


For facial care (page 15)

### ● VI attachment

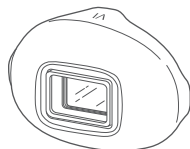

For VI zone care (page 16)

### ● O attachment

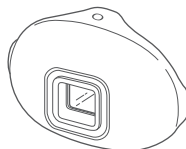

For O zone care (page 16)

- V design tool
- AC adapter
- User's Manual

The design and specifications of the device and attached accessories are subject to change for improvement purposes without notice.

## Functions

### Safety function

#### ● Touch sensor

You can flash your skin once the device detects that the treatment window is fully touching on your skin.

### Auto Power Off

About ten minutes after the device has been powered on, it will be automatically shut off after two beeps. To resume using the device, wait at least ten minutes before proceeding.

### Cooling feature

After the power is turned off, the fan in the device keeps operating for approx. three minutes. It is a function to cool the device.

### Waterproof

This product (excluding the AC adapter) is IPX7 level waterproof\* (JIS waterproof protection class).

- \* The product is not designed to use underwater, so please be careful when using the device in bathrooms or around water.
- Never charge the device when it is wet.
- The AC adapter is not waterproof.

### Charging Guide Display


The device tells you that the battery level is low when the power is turned off. When "L" is shown on the LEVEL indicator and the device beeps six times, it means the battery level is low. Charge the battery before the next use.

## Modes

This product enables you to use the following modes to treat your skin by fitting the supplied attachments to the device (page 22).

### BODY Mode

No attachment

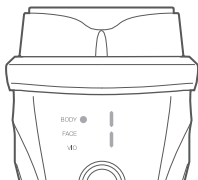

### FACE Mode

Use the facial attachment



<Recommended areas of use>

#### ● Arms



#### ● Legs



#### ● Underarms



#### ● Bikini line



\* Area outside panty line

After setting the facial attachment, the device will automatically switch to FACE Mode.

Be sure to fit the facial attachment to the device and use FACE Mode when applying treatment on your face.

Start using the device from level 1.

<Recommended areas of use>

#### ● Face

(cheeks)



#### ● Face (around the mouth)



\* The device cannot be used on the lips.


#### ● Fingers





### VIO Mode

Use the  
VI or O attachment

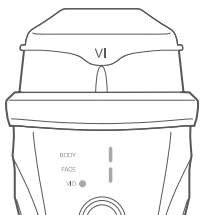


After setting the VI or O attachment,  
the device automatically switches to VIO Mode, flashing the  
light which is appropriate for sensitive skin areas.

When treating the V-, I-, and O-lines, use the attachments  
suitable for each area.

<Recommended areas of use>

#### ● V-line

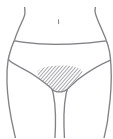


\* Area inside  
panty line

#### ● O-line



Back

#### ● I-line

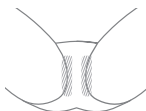


\* Do not flash the light inside  
the I-line, such as on male  
genitalia.

### Areas you should NOT treat

When performing treatment, do not use the device on the following areas.

(Areas shown with  in the figure below)

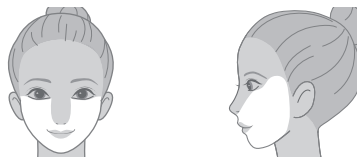

- Anywhere inside the mouth or any mucous membrane
- Eyeballs
- Eyelids (eye sockets)
- Throat
- Bulging moles
- Viral warts
- Areas where you have broken a bone
- Areas of internal bleeding, swelling or scabbing after cosmetic medical procedures
- \* Avoid using the device for 2 weeks after cosmetic medical procedures.  
After this period, consult with a physician and only use the device if the physician confirms there is no problem  
(e.g. hyaluronic acid injections, silicone injections, Botox, thread lifts, double eyelid surgery, HIFU, etc.).
- Dark-colored areas such as scars and spots
- Eyebrows, forehead, and hair (including sideburns)
- ears
- the nape of the neck
- Nipples, areolae
- Inside the navel
- Areas where a large vein or other blood vessels can be seen
- Tattoo or body-painted area

## ■ Before Use

Before carrying out a treatment, be sure to check the following contents.

### Preparing for treatment

- Be sure to remove excess hair at the area to be treated with a handheld shaver or razor.  
A burn injury may result if excess hair is allowed to grow long or remain on the surface of the skin at the area to be treated.  
\* Avoid tweezing or removing hair with wax or other such means.
- Check to make sure there is nothing unusual about the head, treatment window, or attachment. Do not use if damaged or deformed.
  - Be sure to clean the device if it gets dirty. Failure to comply may damage and deteriorate the device or cause burns, injuries, or skin problems.
  - See 'How To Clean' (page 28) for usage instructions.

- Be sure to perform a Flash test before using the device for the first time.


Apply a flash to an area near the intended care area starting with Level 1 in a single flashing. Use BODY Mode for the body, FACE Mode for the face, and VIO Mode for the V-, I-, and O-lines (refer to pages 20-27).

- You can use the device if no abnormality appears on your skin 24 hours after performing the flash test.

- Adolescents must use the device with close attention to safety under the supervision of their parent and, when there is abnormality, stop using the device immediately.


Due to factors related to hormonal balance, damage may be caused to the skin, or the treatment may not yield observable results while the user has not reached full adulthood.

<Treatment test>

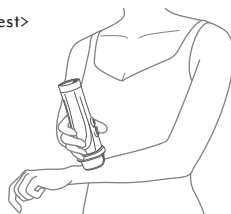

- For hygienic reasons, do not share the V1 and O attachments with your family members or friends.
- Lighting during use
  - Using the device in a very dark environment may cause injury to your eyes.
  - If you want to check the color chart, please use it in bright light.  
(For details on the color chart, please read the included “Safety Precautions”.)
  - Avoid use in the following conditions.
    - An indirectly lit room
    - A very dark room, such as one lit only by miniature bulbs
    - Dark places such as under a table

## Treatment cycle

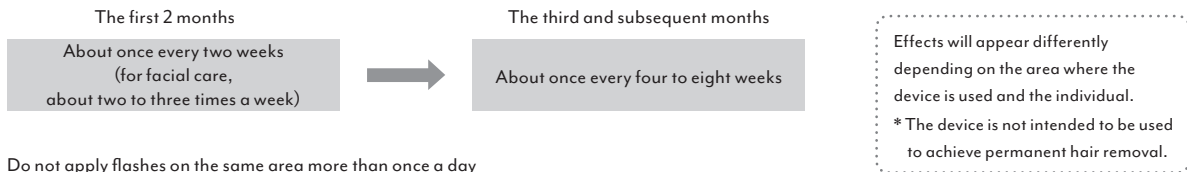


Do not apply flashes on the same area more than once a day  
(As for soft flashing, three high-speed pulses are counted as one flash).

## CAUTION

- Do not engage in excessive use of the device.
- If you wish to use the device in excess of the aforementioned frequency, be sure to use it while checking the condition of the skin. In addition, always adhere to the prescribed frequency of use per day. Failure to comply may cause skin problems.

## Charge the Device

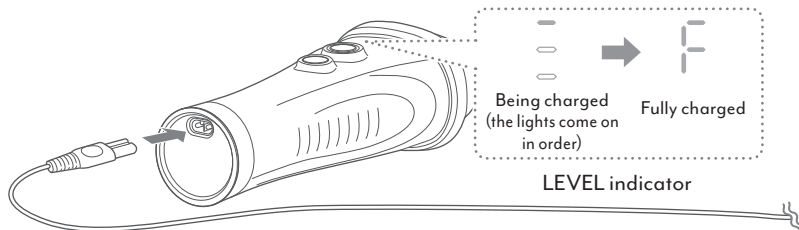

- Charging starts when the AC adapter is connected to the device.
- When the device is being charged, the segments of LEVEL indicator are lit one by one in order. When the device is fully charged, LEVEL indicator displays "F" for 30 seconds and then turns off.
- At an indoor temperature of approx. 25° C, it will take about two and a half hours to reach a full charge, at which the device can operate for about 20 minutes. Actual charging time and operating time may vary depending upon the remaining battery level and the charging environment.
- \* If charging does not begin after approx. 10 seconds, check if the AC adapter connection is loose.

## CAUTION

- **If there is any moisture on or around the main unit or charging terminal, be sure to wipe it dry. Charging while moisture remained could cause an electric shock, short circuit, fire, or a malfunction of the device.**
- **The device cannot be used when it is being charged. Remove the AC adapter before use. If you try to charge the device while the power is on, the power turns off.**
- **Never charge the device in the bathroom.**
- Fully charge the battery before using the device for the first time.
- When using the device for the first time or when the device has not been used for an extended period of time, it is possible that you will be unable to use the device for very long even with a full charge status is displayed. In such cases, recharge the device before using it.
- When the cord is broken, replace the AC adapter with a new one.

## ■ How to Use


### Basic Method of Use

#### 1 Turn on the device

Press and hold the POWER/LEVEL button (approx. 3 seconds). The power turns on after two short beeps.

When MODE indicator and LEVEL indicator are lit and the device emits two short beep sounds, the device is in Standby Mode.

When the device is turned on, the device is set at level 1 of single flashing.

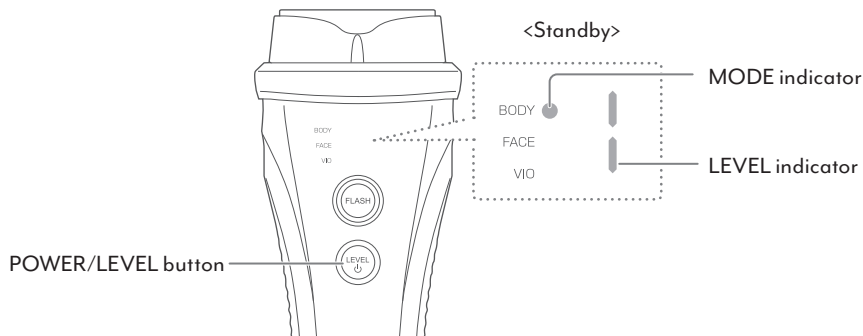




### CAUTION

- Remove the AC adapter before use. The device cannot be used while it is charging.

## 2

## How to set the mode


- Fit an attachment while the device is turned on, the appropriate mode will be selected automatically after one short beep, and the device will start single flashing (page 24) in level 1 (page 22).

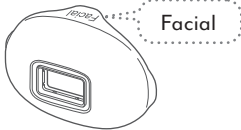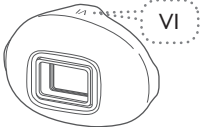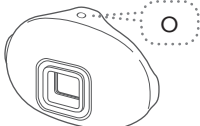Mode	Attachment	MODE indicator
BODY Mode	No attachment	BODY ●      FACE       VIO
FACE Mode	Facial attachment 	BODY       FACE ●      VIO
VIO Mode	VI attachment 	BODY       FACE       VIO ●
	O attachment 	

- Check the mode indicator to confirm the mode.
- Set the right attachment for the area being treated.

## How to connect an attachment

The attachment can only connect to the device in the right direction.

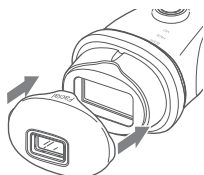


\* The facial attachment is used as an example in the description.

Fit the attachment to the device so that it covers the treatment window, with the attachment's name facing up. Since it is magnetic, it will stay in place if it is covering the treatment window.

- \* Make sure nothing is in between the treatment window and the attachment.
- \* When the attachment is upside down, it cannot be fit to the device.

- If you fit or remove the attachment while the power is on, the mode changes (page 21) after a short beep and the device is automatically set at level 1 (see below) of single flashing (page 24).


## How to remove an attachment

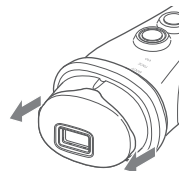

Hold the attachment on both sides and pull it toward you.

## 3 Select a level/irradiation type

Every time you press the POWER/LEVEL button, one beep sounds and the level/irradiation type changes.

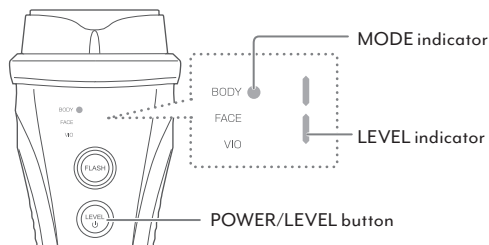

The levels/irradiation type changes as follows.

Every time the irradiation type changes, the device emits a long beep sound followed by a short beep sound.

- Single flashing/continuous flashing (BODY Mode, FACE Mode)



- Single flashing/Soft flashing (VIO Mode)



Check the attached color chart and choose a level based on your skin color.


## CAUTION

- Use only level 1 upon first use and until you are accustomed to use.

## 4

## How to apply treatment

Common to all modes

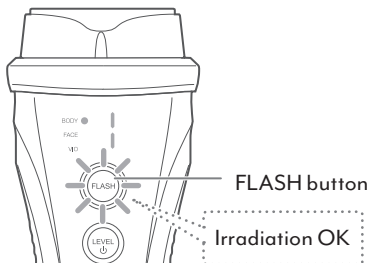

\* The illustration uses BODY Mode as an example.

Put the treatment window on the skin after two short beeps.

When the white light of the FLASH button turns on, press the FLASH button and apply flashes to the skin.

See 'Tips For Your Treatment' (page 25) for tips on performing treatments for different areas.

\* The product is not designed to use underwater, so please be careful when using the device in bathrooms or around water.

\* Wear sunglasses if you feel that the light is too bright during your treatment. Failure to comply may cause damage to your eyes or skin.


• It may take longer to prepare for treatment depending on the level.

• This product is designed to prevent the treatment head from producing excess heat by adjusting the flashing interval, so that the head does not cause damage to users by heating up when continuous flashing is performed continuously.

The flashing interval varies depending on how hot the head is and the environment in which the device is used.

## CAUTION

- Remove the head from the skin after applying a flash to one area and move the head little by little to apply flashes to different areas.
- Since the skin around the face contour is susceptible to getting caught in the treatment window, be careful not to press too strongly while flashing.
- Always wait at least 10 minutes after turning off the power before using the device again. Continual use may cause a malfunction of the device.
- Be careful not to let the filter (transparent glass) touch the skin.

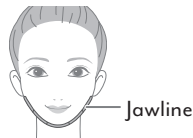



● **Single flashing (once) <can be used in all the modes> MODE indicator is lit**

Press the FLASH button to irradiate a single flashing while shifting the position of the device.

- There are a few seconds of charging time after every single flash (the time varies depending on the level).

● **Continuous flashing <BODY Mode, FACE Mode> MODE indicator flashes**

When you press the LEVEL button when the device is set at level 5 of single flashing, the level changes to level 1 of continuous flashing.

MODE indicator changes to flashing. Choose among levels 1 to 5

(When you press the LEVEL button when the device is set at level 5 of continuous flashing, the level changes to level 1 of single flashing).

- Flashes are emitted while the FLASH button is pressed.
- If you put the treatment window on the skin and press the FLASH button again after continuous flashing is finished, you can perform continuous flashing again (Up to 10 continuous flashes can be emitted in BODY Mode and up to 5 continuous flashes can be emitted in FACE Mode).
- When the treatment window is removed from the skin for 15 seconds, the device beeps and is reset to the state before the first pulse was emitted.




● **Soft flashing <VIO Mode only> MODE indicator flashes**

When you press the LEVEL button when the device is set at level 5 of single flashing, the level changes to level 1 of soft flashing.

MODE indicator changes to flashing. Choose among levels 1 to 5

(When you press the LEVEL button when the device is set at level 5 of soft flashing, the level changes to level 1 of single flashing).

- Press the Flash button to flash light while moving the device over your skin.
- Three high-speed soft continuous flashes are applied to the same area.
- It is recommended for those who do not like the instant stimulation caused by single flashing and those who are sensitive to heat.

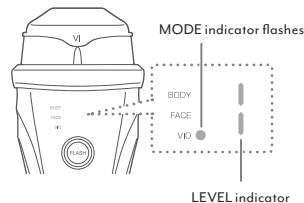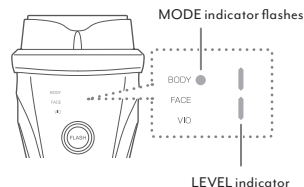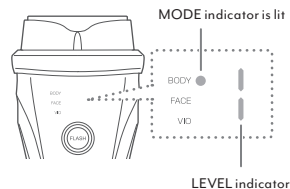

## Tips For Your Treatment

Keep the following points in mind when performing treatment.

### ● Arms/legs

Hold the treatment window perpendicular to the skin.

### ● Underarms

Perform treatment while stretching out the skin with your fingers.

• Before performing treatment, carefully wipe away any antiperspirant or other cosmetic products on the skin.

### ● Bikini line

Use the device at a low level on areas where there is pigmentation or areas with darker skin tone.

\* Be sure to start from level 1 and adjust as you become more accustomed to the device.

### ● Face

Refrain from applying flash treatments without first removing foundation and other makeup products that have been applied to the skin.


Failure to comply may cause a burn injury or skin problems.

• Use the device with a degree of intensity sufficient for you to perceive a slight warmth.

• When using the device on a narrow area such as mustache or beard, or a area with a sharp angle, perform the treatment while stretching out the skin with your fingers.

• Refrain from tanning while undergoing treatment and avoid direct sunlight after treatment.

• When applying flashes, press the treatment window firmly against the skin and be careful not to expose your eyes to the flashes.

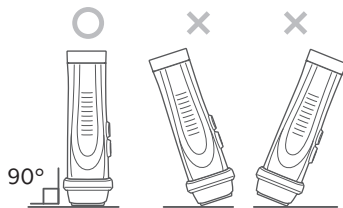


## ● VIO area

### V-line


Perform treatment while stretching out the skin with your fingers.

When you style your pubic hair using the attached V designing tool, shave the area inside the V-line completely and apply flashes.

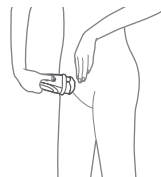


<About the V designing tool>


Triangle

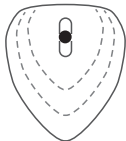


Small triangle

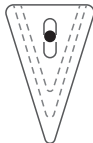

Oval

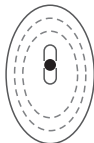

Square

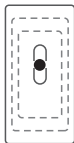

### ● Tab

You can choose among the four designs and style your pubic hair inside the V-line based on it. Each design is adjustable to three sizes.

<How to Use>

1 Choose the shape and decide on the size.

2 Cut the tool to the right size for you along the dotted line with scissors.

\* Adjust the size from large to small

3 Shave the area inside the V-line.

\* Even if you want to keep your hair, make sure to shave the area completely.

4 Push the tab to one side.

5 Press the tool against the V-line area with the hand that is not used to hold the device.

6 Use the device to apply flashes to the areas around the tool.

\* From the second treatment onward, perform the treatment by pressing the same tool you used for the first treatment to the same area.

## I- and O-lines

<Recommended usage>

Put a mirror on the floor and get down on one knee over the mirror as shown in the figure.

Check the area to apply flashes to in the mirror.

- \* Check whether the area to treat is completely shaved.
- \* After checking the area in the mirror, apply flashes without directly looking into the mirror.
- \* Do not apply flashes directly to the mucosa.
- \* Do not flash the light inside the I-line, such as on male genitalia.



Mirror

## 5 Turn off the device

- Press and hold the POWER/LEVEL button for approx. 3 seconds.
- Please confirm that the fan in the device has stopped (approx. three minutes).
- After each use, refer to 'After Treatment' (page 28) and remove the remaining sebum, serum and other particles that have adhered to the device.

If your skin remains hot feeling after use, cool your skin down sufficiently by applying a cold towel.

Your skin may feel hot after carrying out a treatment depending on your physical constitution and method of use. If you are concerned about the hot flashes on your skin, cool it down by applying a ice pack or cold towel to the area. Consult a physician if the abnormality does not subside.

## CAUTION

- Be sure to use the mode that is appropriate to the area being treated. Failure to comply may cause a burn injury.
- Ensure that neither you nor anyone else looks directly into the treatment window of the device. Also ensure that no one nearby looks at the treatment window. Failure to comply may cause injury to your eyes.
- Moisturize your skin daily to prevent your skin from becoming dry. If your skin is dry, the device might not work properly.
- When treating the face, press the treatment window firmly against the skin. If the treatment window is not pressed firmly against the skin, the flashing light may become overly bright.

## After Treatment

### How To Clean

Please carefully clean the parts of the device that come into contact with the skin every time you use it and keep the device clean at all times.

Failure to comply may cause skin problems.

Regularly clean the device also to prevent it from failing or malfunctioning.

Be sure to turn the power off and wipe off moisture after use.

#### Flash irradiation window

- 1 Use a soft piece of cloth or cotton swab to remove stains and bits of loose hair.

\*Wipe away fingerprints and other stains. Failure to comply may cause the filter to become scorched, a diminishment of flash performance, a malfunction of the device, or skin problems.

#### Main unit

- 1 Wipe away any surface dirt with a dampened, well wrung cloth.

- For stubborn stains, wipe with a sponge or cloth soaked in a diluted neutral detergent.
- \* If there is any moisture on or around the charging terminal, be sure to wipe it dry.

#### Attachment

- 1 Remove an attachment from the device and wash it under running water.
- 2 Wipe off moisture with a soft cloth.

#### V design tool

- 1 Wash the tool under running water.
- 2 Wipe off moisture with a soft cloth or tissue.

Do not wipe the device with solvents such as alcohol, thinner, benzene, nail polish remover, kerosene, chemicals, or polishing powder.  
Failure to comply may cause damage to the surface of the device.



## Storing the Device

- Always remove excess water and dry thoroughly before storing. Leaving the device wet may cause a malfunction or corrosion.
- If not using the device for a long period of time, clean the device before storing it.
- Turn the device off when not in use.
- Avoid using or storing the device near a heater, in a high-temperature/high-humidity location, locations where there is exposure to direct sunlight, or unstable locations.
- Keep the device out of the reach of children.
- Do not keep the device in wet locations such as the bathroom.

## Maintenance and Inspections

- If you have not used the device for an extended period, inspect it carefully and check that it works properly and safely before using it again.
- If the device is malfunctioning, go to 'Troubleshooting' (page 30) and attempt to identify the problem before arranging to have the device repaired. If the problem persists, contact your local distributor, making sure that you have your purchase details ready when you call.
- If the device malfunctions, appears damaged, or needs service, please contact your local distributor.

## Additional Information

### Troubleshooting

If you notice an abnormality during use, check for the following symptoms.

Symptom	Possible causes and solutions
The device does not turn on.	● The device has run out of power. Charge the battery (page 19).
	● Press and hold the POWER/LEVEL button for approx. 3 seconds to turn on the device.
	● The device cannot be turned on while being charged. Remove the AC adapter.
The device cannot be operated.	● Read the “Error display” (page 32).
The device won't stay on.	● The device is set to shut off automatically (page 14). To resume using the device, wait at least ten minutes before proceeding.
	● The battery is running low. Charge the battery (page 19).
No flash pulse is emitted.	● Preparations for flash emission have not been completed. Wait until 2 short beeps sound.
	● Ensure that the treatment window is fully touching over your skin and wait for the lamp of FLASH button to turn on before removing your hair (page 23).
	● The glass tube is at the end of its life. Read the “Error display” (page 32).
No flash is emitted (single flashing).	● Press the FLASH button (page 23).
	● Press the treatment window firmly on your skin.

**Symptom****Possible causes and solutions**

**No flash is emitted  
(continuous flashing).**

- Press the treatment window firmly on your skin.
- Flashes are emitted while the FLASH button is pressed (page 24).
- Check whether MODE indicator is flashing (page 22).

**The flash pulse feels  
hot or painful.**

- Press the POWER/LEVEL button to lower the level (page 22).
- There is excessive hair. Make sure to remove excess hair with a handheld shaver or razor.
- Review the areas on which use of the device is prohibited (pages 5 and 16).
- Your skin has entered inside the treatment window. Do not press your skin strongly.
- Do not apply flashes to the same area more than once (As for soft flashing, three high-speed continuous flashes are counted as one flash).

**Flash pulses are bright.**

- The treatment window is not in full contact with the skin. Have the treatment window come into firm contact with the skin.
- Please wear sunglasses.

**There is a buzzing sound.**

- It is a sound generated by the fan operating in the device. This is not an abnormality or malfunction.

**The flash preparation time  
is long.**

- The higher the set level is, the longer the preparation time for per flash.
- When the treatment head heats up, the preparation time per flash becomes longer.



## Error display

When an error occurs, error information is displayed on LEVEL indicator of the device.

Check the error in question as outlined below and eliminate the cause of the error accordingly.

Error indication	Solutions
E and H light up alternately	The power is turned off automatically after five short beeps. The head is hot, allow the device to rest for about 30 minutes.
E and I light up alternately	The power is turned off automatically after five short beeps. As the glass tube is at the end of its life, contact your local distributor.
Either E or 2 to 5 light up alternately	The power is turned off automatically after five short beeps. Inquire with your local distributor (page 36).

## Frequently Asked Questions

Q

The flash pulses feel hot and redness or pain is recognized. Is this normal?

A

- Start using the device from level 1.
- Once you become comfortable with this level, gradually raise the level. Perceiving a slight measure of heat from flash pulses is normal and does not cause any skin problem.
- This problem can be attributed to the presence of long bits of excess hair or to the treatment of suntanned skin with flash pulses.
- Use the device after getting rid of excess hair with a handheld shaver or razor and on skin that has not been suntanned.
- If redness or pain is manifested, discontinue using the device immediately and cool down and moisturize the skin. If there is no amelioration, discontinue using the device and see a physician.  
\* Be sure to use the facial attachment when treating the face and the VI attachment or O attachment when treating the private parts. Failure to comply may cause a burn injury.
- Be careful not to let the filter (transparent glass) touch the skin. When the device is being used, the filter becomes extremely hot and may cause skin burns.
- Do not press the device too firmly against the areas where the skin is soft or there is a lot of fat under the skin. Be careful when using the device to treat the rounded parts, such as elbows and knees.

Q

How often can I use the device per day?

A

Do not flash-illuminate the same area more than once a day. Intense heat-generated stimulation may give rise to skin problems.

Q

Immediately after a flash pulse is discharged, I can smell a burning odor. Is this normal?

A

The burning odor comes from the heating of excess hair caused by the flash pulse. If a burning odor is emitted by the device itself, immediately discontinue use and turn the device off.

Q

Can I use the device in the bathroom?

A

Though the device can be used in the bathroom, it cannot be used underwater. In addition, do not charge the device in the bathroom.

Q

Can I share the device with my family members or friends?

A

The device can be shared for BODY Mode that does not require any attachment or treatment with the facial attachment. When sharing the device, be sure to clean the device and attachment before use. Do not share the VI attachment and O attachment with anyone.

Q

Can I use the device on the skin on which the cooling/moisturizing gel is used?

A

Do not use the device because it may cause malfunction. Use the gel when there is redness or pain after use.

## Product Specifications

<b>Rated power supply</b>	DC9V 2.0A
<b>Power Consumption</b>	Approx. 8.0 W (when charging)
<b>Outer Dimensions</b>	Approx. 82mm × 65mm × 203mm (W×D×H)
<b>Weight</b>	Approx. 450 g (main unit only)
<b>Water Resistance Class</b>	IPX7
<b>Materials</b>	Device: Polycarbonate resin, nylon resin, ABS resin, glass Facial attachment: Polycarbonate resin, stainless steel (treatment window), glass, silicone VI attachment: Polycarbonate resin, stainless steel (treatment window), glass, silicone O attachment: Polycarbonate resin, stainless steel (treatment window), glass, silicone
<b>Country of Origin</b>	China
<b>Accessories</b>	One facial attachment One VI attachment One O attachment One V designing tool One AC adapter User's Manual

## ■ After-sales service

### Rechargeable battery

- The rechargeable battery is a consumable item (It is not covered by warranty).
- This appliance contains batteries that are non-replaceable.
- If the device no longer operates for more than a short period of time even after the battery has been recharged, the battery may have reached the end of its service life (the battery is expected to last 300 charge cycles).
- For battery replacement, please contact your local distributor.

### DISPOSAL INFORMATION

- Follow your local government rules on separation of trash for proper disposal of this product.
- The battery must be removed from the appliance before it is disposed of.
- The lithium-ion battery used in this product is a valuable recyclable resource.



Li-ion

### Contacts for customers outside Japan

- Please contact your local distributor or visit YA-MAN Global websites (<https://www.ya-man.com/en/>).

## ■ Marking Description



This marking indicates that this product should not be disposed of with other household wastes throughout the UK and EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources.

To return your used device, please use the return and collection systems or contact the retailer where the product was purchased, they can take this product for environmental safe recycling.




This marking indicates that it is a detachable supply unit.  
This product must be charged only using the included detachable power supply.


## How to remove the built-in battery

- Be sure to turn the power off before removing the battery.
- Use a screwdriver to remove the battery. Be careful not to injure yourself.
- When removing the battery, be careful not to damage it with a hard object.\* Failure to comply may cause the generation of heat or a fire.
- Follow your local government rules for proper disposal or recycling of this product.


**1** Remove the ring part fixed at the base of the head using an object with a sharp tip.

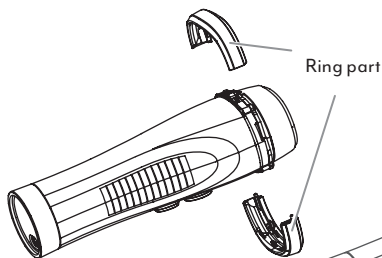

**2** Remove the four screws and take the internal unit out of the case.

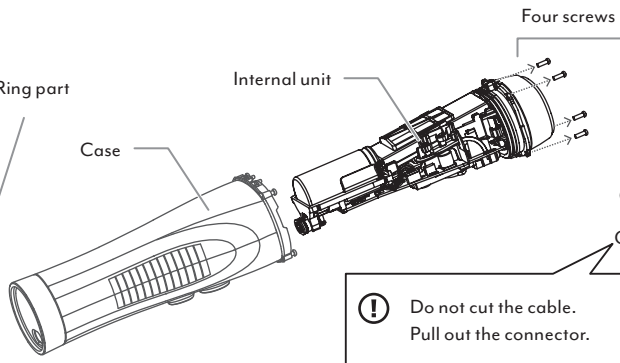

**3** Remove the four screws and then remove the battery.

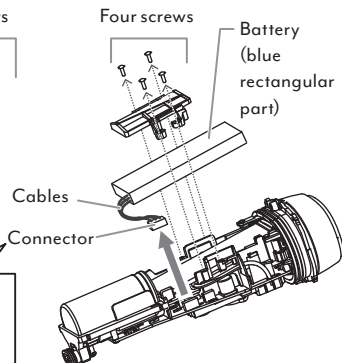

Above figures are intended to be used for disposing of the product, not repairing it.

**CAUTION:** Disassembling the device by yourself may cause a failure.

**Memo**

